
BiteCalm

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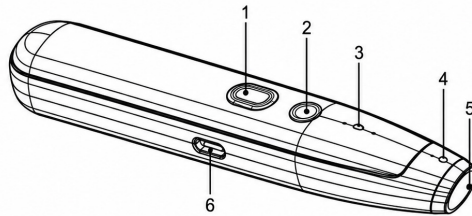BiteCalm

Electronic Insect Bite Relief Device - SB-1316

4 OPERATING MODES

Choose the mode that best suits your needs and skin sensitivity.

FAST RELIEF

Helps relieve itching and discomfort.

EASY TO USE

Intuitive controls with an easy-to-read display.

COMPACT DESIGN

Lightweight and convenient to carry anywhere.

The Cangaroo Compact Electronic Insect Bite Relief Device is a practical solution for quickly soothing the skin after insect bites and stings – without the need for creams or chemicals.

Featuring a modern design, easy operation, and 4 operating modes, the device allows you to choose the most suitable intensity according to your skin sensitivity and personal comfort. It offers precise temperature control with 4 temperature settings and adjustable treatment duration, providing fast relief from itching, pain, and swelling caused by insect bites and stings. Suitable for adults and children over 6 years of age.

PRINCIPLE OF ACTION

When an insect bite or sting occurs, the insect's saliva triggers a natural defense response in the body. Immune cells are activated, histamine is released, and signals are sent to nerve receptors in the skin.

This immune response may cause redness, swelling, irritation, and itching. **BiteCalm** helps soothe the skin reaction through a method known as local thermotherapy, which is based on two key principles:

1. Heat as a Counter-Stimulus

Heat helps inhibit the transmission of itch signals to the nerve cells in the skin, thereby relieving the sensation of itching.

2. Heat Reduces the Effects of Histamine

Studies suggest that heat can stimulate beneficial physiological responses, including a reduction in histamine activity, which may help alleviate itching and skin irritation.

HOW IT WORKS

Temperature Selection

The device offers 4 temperature settings:

- Level 1: 40°C
- Level 2: 43°C
- Level 3: 46°C
- Level 4: 51°C

Press and hold the temperature selection button for 2 seconds to choose the desired temperature setting. Once a temperature setting has been selected, the device will automatically switch off after 60 seconds if no treatment is started. Press the treatment button to start the treatment. The device will automatically switch off after 30 seconds when the treatment cycle is completed.

CHARGING

Connect the USB Type-C cable to a power adapter, computer, or other USB power source.

- The temperature indicator flashes while charging.
- When the battery is fully charged, the indicator remains steadily lit.

SAFETY WARNINGS

- Use the device only for its intended purpose – relieving symptoms

caused by insect bites and stings.

- Do not use on broken, irritated, or inflamed skin.
- Do not apply to sensitive areas around the eyes or on mucous membranes.
- The device heats up during use; use only as directed.
- Do not use if you have impaired sensitivity to heat or are unable to perceive temperature correctly.
- Children should use the device only under adult supervision.
- Discontinue use if you experience severe discomfort, pain, or skin irritation.
- Do not immerse the device in water or use it with wet hands.
- Keep out of reach of children.
- In the event of a severe allergic reaction or significant swelling, seek immediate medical attention.

CLEANING

- Clean the contact surface with a soft, slightly damp cloth after each use.
- Do not immerse the device in water.
- Ensure the device is completely dry before storage.

Read all instructions and warnings carefully before using this product. KEEP FOR FUTURE REFERENCES.

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Manufacturer and Importer: Moni Trade Ltd.