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EN

INSTRUCTIONS FOR USE

IMPORTANT! READ CAREFULLY AND KEEP FOR FUTURE REFERENCE!

TO ASSURE SAFE USAGE OF THIS PRODUCT FOLLOW ALL WARNINGS, INSTRUCTIONS AND RECOMMENDATIONS GIVEN IN THIS USER MANUAL!

1. INTENDED USE

The Kangaroo electronic scale is designed to meet the needs of the entire family. The combined function of the scale allows: 1. Precise measurement of weight of babies in lying position with the help of removable tray. 2. User-friendly weight measurement of adults and older children in standing position. The device has been tested and complies with the requirements of European standards EN IEC 55014-1:2021 and EN IEC 55014-2:2021.

The scale has the following features:

- Large, easy-to-read LCD display.
- Automatic and manual hold function for recording the weight of wriggling babies.
- Tare weighing.
- 180 kg weight capacity.
- Automatic switch-off.
- Adjustable between weight units (kg/lb/oz/lb:oz).

PACKAGE CONTAINS:

- Electronic Scale - 1 pc.
- User Manual - 1 pc.
- 3 x 1.5V AAA batteries - **NOT INCLUDED IN THE PACKAGE.**

2. RECOMMENDATIONS AND WARNINGS FOR SAFE USE OF THE PRODUCT

WARNING!

- Use the appliance only for its intended purpose and in accordance with the instructions and recommendations described in this manual. Using the appliance in a way other than intended may result in damage and/or injuries. Any damage caused by failure to follow these instructions will void the warranty.
- This product is intended for household use **ONLY** for weighing infants and small children. It is not intended for commercial or industrial use.
- The surfaces of this scale are not intended to come into contact with food.
- This device is a precision measuring instrument. Please handle it with care. Do not drop the scale itself or place other objects on it that could cause impact.
- Protect the scale from contact with sharp objects, moisture, dust, chemicals, toiletries, liquid cosmetics, significant temperature fluctuations, and sources of heat.
- Repairs or servicing of the device must be carried out only by authorized technical personnel.
- Do not press the buttons forcefully or use sharp objects to operate them.
- Do not place objects on the scale when it is not in use.
- Strong electromagnetic fields (e.g., from a mobile phone) may affect the accuracy of the scale.
- Remove the transport safety lock before using the scale.
- The maximum capacity of the scale is 180 kg (396 lbs).
- The factory-set unit of measurement for this scale is kilograms.
- Always use the scale on a horizontal, solid, and even surface.
- **WARNING!** Never leave a child unattended on the scale.

- Keep the scale’s packaging out of the reach of children – risk of suffocation!

3. IMPORTANT BATTERY GUIDELINES

WARNING!

- Keep batteries out of the reach of children due to the risk of accidental ingestion.
- Always insert the batteries correctly, observing the polarity (+/-) markings inside the battery compartment.
- When replacing batteries, always replace all of them at the same time. Do not mix old batteries with new ones.
- Do not dispose of depleted batteries with household waste. Dispose of them at designated collection points. Batteries are recyclable.
- If a battery leaks, wear protective gloves and clean the battery compartment with a dry cloth. If battery fluid comes into contact with your skin or eyes, rinse immediately with water and seek medical attention.
- If the device will not be used for an extended period, remove the batteries.
- Protect batteries from excessive heat and never dispose of them in fire.
- Use only AAA-size batteries. Do not use rechargeable batteries.
- If battery fluid comes into contact with your skin or eyes, rinse thoroughly with water and seek medical assistance.

4. TECHNICAL SPECIFICATIONS

Scale batteries: 3x 1.5V AAA (not included in the package); **Measuring range:** 0.5 kg - 180 kg (1.1 lb – 396 lb); **Measuring accuracy:** ±0.1 kg (±0.2 lbs); **Resolution:** 0.01 lb (5 g) increments (0 to 44 lb / 0 to 20 kg); 0.1 lb (50 g) increments (44 to 396 lb / 20 to 180 kg); **TEMPERATURE RANGES:** **Operating temperature:** 5°C to 40°C (41°F to 95°F), **relative humidity:** 20 - 85%, **storage temperature:** 20°C to 60°C (68°F to 140°F), **relative humidity:** 20 - 90%; **Dimensions with tray:** 515 x 327 x 71 mm (20.3 x 12.9 x 2.8 inches); **Dimensions without tray:** 302 x 302 x 27 mm (11.9 x 11.9 x 1.1 inches); **Net weight of the scale (with tray and batteries):** 2.01 kg (4.4 lb).


5. DEVICE DESCRIPTION

FRONT SIDE - FIGURE 1

1. LCD/Hold button.
2. Display.
3. POWER/TARE button.
4. Movable tray.

BACK SIDE - FIGURE 2

5. “Weight Selector” button (Unit selection).
6. Battery compartment.
7. Fixators for securing the tray (tray tabs).

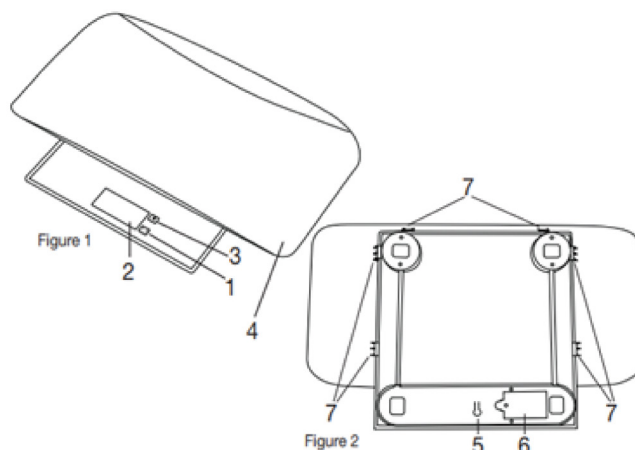REMARK: SELECTING THE UNIT OF WEIGHT (WEIGHT SELECTOR BUTTON): When supplied to the customer, the baby scale is set to weigh in kilograms “kg”. You can find the unit button (weight selector) on the rear of the baby scale. You can also use this to change the unit of weight to pounds “lb” or ounces “oz”.


6. INSERTING THE BATTERIES

Before using the baby scale, you should first insert three (3) 1.5 V AAA batteries into the battery compartment. To do this, open the compartment cover on the bottom of the scale, insert the batteries while observing the correct polarity. Remove the battery compartment cover, located under the scale, by unscrewing the fixing screw counterclockwise and lift the cover away by the tab fixator. Then close the cover carefully. Make sure that the batteries are inserted the right way round (+/-), as shown in the battery compartment.

NOTE: The scale is equipped with a low-battery indicator. When the symbol “Lo” appears on the display, the scale will switch off automatically. The batteries must then be replaced with new ones.

7. HOW TO INSTALL THE TRAY

Slide the baby tray onto the scale by aligning the slots on the underside of the tray with the sides of the scale and pushing the tray onto the scale until the tabs on the bottom of the tray contact the bottom edge of the scale.

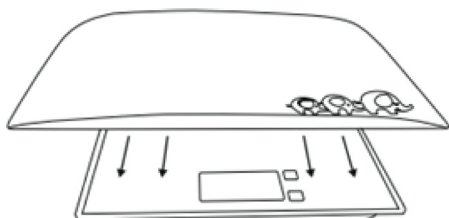8. OPERATION OF THE DEVICE

8.1 WEIGHING

Place the scale on a horizontal, solid, flat, and vibration-free surface. Placing the scale on a soft and/or uneven surface (e.g., carpet, linoleum, etc.) may result in inaccurate readings.

- Press the **POWER/TARE** to switch on the scale. The scale will perform a self-diagnostic check. The scale is ready for weighing as soon as the display shows “0.000 kg.”
- Press the **LCD/HOLD** button to set the desired display orientation.
- Carefully place the baby in the center of the weighing platform (scale tray). The weight is displayed on the screen and then automatically turns off.
- Press the **POWER/TARE** button again to perform another measurement.
- To switch off the scale manually, press and hold the **POWER/TARE** button for 3 seconds.

8.2 HOLD FUNCTION (WEIGHT READING HOLD)

The scale is equipped with a built-in hold function (displaying the average reading). This feature allows accurate weighing of the baby, even when the baby is moving and not lying calm on the weighing platform.

The hold function (“HOLD”) offers two modes: Automatic Hold (AUTO HOLD) and Manual Hold (MECHANICAL HOLD)

AUTOMATIC HOLD:

- Switch on the scale. The scale is ready for weighing as soon as the display shows “0.000 kg.”
- Carefully place the baby in the center of the weighing platform (scale tray).
- The AUTO HOLD function (activated by default) determines the baby’s weight as soon as the child remains sufficiently still to allow an accurate measurement. It is important not to touch the child during this process, as doing so may affect the accuracy of the measurement.
- As soon as the baby’s weight is determined, the “HOLD” symbol stops flashing, and the weight is displayed on the screen for 2 minutes.
- The measured weight remains stored in the scale, and the device will automatically switch off after 2 minutes. To turn it on again, press and hold the **POWER/TARE** button for 3 seconds.

MANUAL HOLD (MECHANICAL HOLD):

If the baby is excessively restless or moving continuously, you can activate the manual hold function by pressing the **LCD/HOLD** button. After 5 seconds, the “HOLD” symbol will stop flashing, and the average weight will be displayed on the screen for 2 minutes.

8.3 TARE FUNCTION

The scale features an automatic TARE function for weights ranging from 0 to 500 grams. This ensures that during subsequent measurements, only the actual weight of the baby will be displayed. The weight will remain visible for 5 seconds before resetting to zero. If you place anything over 500 grams on the scale, briefly press the **POWER/TARE** button to execute the Tare function manually. “-----” will briefly appear on the LCD Display and the weight resets to zero.

8.4 MEASURING THE WEIGHT OF ADULTS AND OLDER CHILDREN IN STANDING POSITION

- Slide the tray off the scale (if attached).
- Press the **POWER/TARE** to switch on the scale.
- Press the **LCD/HOLD** button to orient the LCD Display so it is readable while you are standing on it with the LCD Display toward the front.

NOTE: Power/Tare Button will become inoperable when the scale is being used in this manner to prevent inadvertent activation from a user’s feet if they are standing on top of the LCD/Hold Button.

- The scale will perform a self-diagnostic check. The scale is ready for weighing as soon as the display shows “0.000 kg.”
- Stand motionless on the scale with your feet apart until the weight stabilizes and appears in the LCD Display.
- The scale will shut off by itself after about 30 seconds.
- To reactivate the scale, tap your foot on the scale.

NOTE: To reactivate the **POWER/TARE** Button, reorient the LCD Display by briefly pressing the **LCD/HOLD** button.

9. CLEANING AND MAINTENANCE

- The accuracy of the measurements and service life of the device depend on its careful handling: The baby scale should be cleaned from time to time.
- This should be done using a damp cloth and, if necessary, a small amount of detergent.
- Do not submerge the baby scale in water. Do not rinse the baby scale under running water. Do not clean the baby scale in the dish washer.
- Never use harsh solvents or cleaning products!
- If you do not plan to use the baby scale for a long period of time, we recommend that you store it in the original packaging in a dry environment.
- During storage, do not place any objects on the baby scale.

- Store the baby scale away from children and pets.
- Remove the batteries from the baby scale.

10.TROUBLESHOOTING GUIDE

DISPLAY	CAUSE	SOLUTION
EEEE	The surface on which you have placed the scale is not stable and/or flat.	Place the scale on a flat stable surface without carpet. Retare the scale and repeat the measurement.
	The maximum scale weight of 180 kg (396 lbs) has been exceeded.	Do not exceed the limit on the scale.
Incorrect weight is displayed	The surface on which you have placed the scale is not stable and/or flat.	Place the scale on a flat stable surface without carpet. Retare the scale and repeat the measurement.
	The baby kicks its legs and moves too much.	Try to calm the baby so you can get accurate measurements.
	Incorrect scale reset. Incorrect zero point of scale.	Turn the scale off, turn it back on and wait for "0.000 kg" to appear on the display. Repeat the measurement.
Lo	The batteries are empty.	Replace the batteries in the baby scale.

WARNING!

If the device does not operate even after performing the above procedures, please contact the store where you purchased it or the nearest authorized distributor. Never turn on the device if any part of it is damaged or not functioning properly. The device should only be repaired by an authorized customer service representative or supplier.

By law, empty batteries must be disposed of at designated locations, recycling points or electronics retailers. The codes below are printed on batteries containing harmful substances:

- Pb = Battery contains lead
- Cd = Battery contains cadmium
- Hg = Battery contains mercury

For environmental reasons, do not dispose of the device in household waste at the end of its service life. Dispose of the device at a collection or recycling point.

This appliance complies with the following documents:
 Directive 1999/5/EC on radio equipment and telecommunications terminal equipment.
 Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment.

Manufactured for CANGAROO in PRC;
Manufacturer and Importer: Moni Trade LTD
Address: Bulgaria, Sofia, Trebich, Dolo str. 1 Tel.: +359 2 936 07 90 Website:
www.cangaroo-bg.com

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