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EN

INSTRUCTIONS FOR USE

IMPORTANT! READ CAREFULLY AND KEEP FOR FUTURE REFERENCE!

TO ASSURE SAFE USAGE OF THIS PRODUCT FOLLOW ALL WARNINGS, INSTRUCTIONS AND RECOMMENDATIONS GIVEN IN THIS USER MANUAL!

1. INTENDED USE

The Kangaroo electronic baby scale helps you track your baby's healthy growth with ease. It is designed to assist in monitoring the healthy development of infants by providing accurate measurements of two essential parameters: the baby's weight and height.. The device has been tested and complies with the requirements of European standards EN IEC 55014-1:2021 and EN IEC 55014-2:2021.

The baby scale has the following features:

- Curved weighing surface.
- Large, easy-to-read LCD display.
- Automatic and manual hold function for recording the weight of wriggling babies.
- Integrated tape measure (150 cm).
- Tare weighing.
- 20 kg weight capacity, 5 g graduation.
- Automatic switch-off.
- Adjustable between weight units (kg/lb/oz/lb:oz).

PACKAGE CONTAINS:

- Baby Scale with Height Measuring - 1 pc.
- User Manual - 1 pc.
- 3 x 1.5V AAA batteries - **NOT INCLUDED IN THE PACKAGE.**

2. RECOMMENDATIONS AND WARNINGS FOR SAFE USE OF THE PRODUCT

WARNING!

- Use the appliance only for its intended purpose and in accordance with the instructions and recommendations described in this manual. Using the appliance in a way other than intended may result in damage and/or injuries. Any damage caused by failure to follow these instructions will void the warranty.
- This product is intended for household use ONLY for weighing infants and small children. It is not intended for commercial or industrial use.
- The surfaces of this scale are not intended to come into contact with food.
- This device is a precision measuring instrument. Please handle it with care. Do not drop the scale itself or place other objects on it that could cause impact.
- Protect the scale from contact with sharp objects, moisture, dust, chemicals, toiletries, liquid cosmetics, significant temperature fluctuations, and sources of heat.
- Repairs or servicing of the device must be carried out only by authorized technical personnel.
- Do not press the buttons forcefully or use sharp objects to operate them.
- Do not place objects on the scale when it is not in use.
- Strong electromagnetic fields (e.g., from a mobile phone) may affect the accuracy of the scale.
- Remove the transport safety lock before using the scale.
- The maximum capacity of the scale is 20 kg (44 lbs / 704 oz), and the measurement increment is 5 g (0.02 lbs / 0.20 oz).
- The factory-set unit of measurement for this scale is kilograms.

- Always use the scale on a horizontal, solid, and even surface.
- **WARNING!** Never leave a child unattended on the scale.
- Keep the scale’s packaging out of the reach of children – risk of suffocation!

3. IMPORTANT BATTERY GUIDELINES

WARNING!


- Keep batteries out of the reach of children due to the risk of accidental ingestion.
- Always insert the batteries correctly, observing the polarity (+/-) markings inside the battery compartment.
- When replacing batteries, always replace all of them at the same time. Do not mix old batteries with new ones.
- Do not dispose of depleted batteries with household waste. Dispose of them at designated collection points. Batteries are recyclable.
- If a battery leaks, wear protective gloves and clean the battery compartment with a dry cloth. If battery fluid comes into contact with your skin or eyes, rinse immediately with water and seek medical attention.
- If the device will not be used for an extended period, remove the batteries.
- Protect batteries from excessive heat and never dispose of them in fire.
- Use only AAA-size batteries. Do not use rechargeable batteries.
- If battery fluid comes into contact with your skin or eyes, rinse thoroughly with water and seek medical assistance.

4. TECHNICAL SPECIFICATIONS

Model: EB524; **Dimensions:** 560 × 330 × 47.5 mm; **Weight:** 1.35 kg (without batteries); **Maximum Voltage:** 4.5 V; **Operating Temperature:** up to 40 °C; **Rated Current:** 15 mA; **Batteries:** 3 × 1.5 V AAA; **Measuring Range:** 0.050 kg – 20 kg (0.11 lb – 44 lb, 1.76 oz – 705 oz).

5. DEVICE DESCRIPTION

1. Display.
2. ON/OFF  button.
3. TARE function button.
4. HOLD function  button.
5. Weighing surface.
6. Measuring tape for height measurement, retractable from the scale.
7. Unit selection button (to change the unit of measurement).
8. Battery compartment with cover.
9. Label with essential technical data.

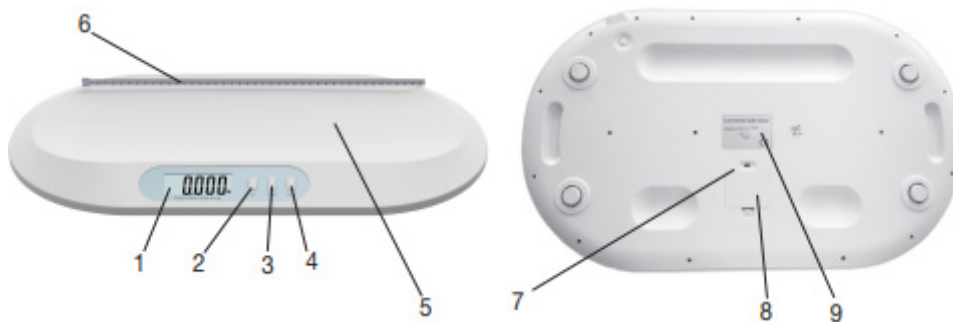
DISPLAY INDICATIONS

1. “TARE” – Indicates that the TARE function is active
2. “HOLD” – Indicates that the HOLD function is active
3. Weight display – Shows the measured weight
4. Unit of measurement – Pound (lb), Ounce (oz), Kilogram (kg).

REMARK: SELECTING THE UNIT OF WEIGHT: When supplied to the customer, the baby scale is set to weigh in kilograms “kg”. You can find the unit adjuster on the rear of the baby scale. You can also use this to change the unit of weight to pounds “lb” or ounces “oz”.

5. INSERTING THE BATTERIES

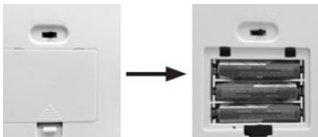Before using the baby scale, you should first insert three (3) 1.5 V AAA batteries into the battery compartment. To do this, open the compartment cover on the bottom of the scale, insert the batteries while observing the correct polarity, and then close the cover carefully. Make sure that the batteries are inserted the right way round (+/-), as shown in the battery compartment.

NOTE: The scale is equipped with a low-battery indicator. When the symbol “Lo” appears on the display, the scale will switch off automatically. The batteries must then be replaced with new ones.

7. OPERATION OF THE DEVICE

7.1 WEIGHING

Place the scale on a horizontal, solid, flat, and vibration-free surface. Placing the scale on a soft and/or uneven surface (e.g., carpet, linoleum, etc.) may result in inaccurate readings.

- Press the ON/OFF button to switch on the scale. The scale will perform a self-diagnostic check. The scale is ready for weighing as soon as the display shows “0.000 kg.”
- Carefully place the baby in the center of the weighing platform. The weight will be displayed on the screen for 2 minutes, after which the scale will automatically switch off.
- Press the ON/OFF button again to perform another measurement.
- To switch off the scale manually, press and hold the ON/OFF button for 3 seconds.

7.2 HOLD FUNCTION (WEIGHT READING HOLD)

The scale is equipped with a built-in hold function (displaying the average reading). This feature allows accurate weighing of the baby, even when the baby is moving and not lying calm on the weighing platform.

The hold function (“HOLD”) offers two modes: Automatic Hold (AUTO HOLD) and Manual Hold (MECHANICAL HOLD)

AUTOMATIC HOLD:

- Switch on the scale. The scale is ready for weighing as soon as the display shows “0.000 kg.”
- Carefully place the baby in the center of the weighing platform.
- The AUTO HOLD function (activated by default) determines the baby’s weight as soon as the child remains sufficiently still to allow an accurate measurement. It is important not to touch the child during this process, as doing so may affect the accuracy of the measurement.
- As soon as the baby’s weight is determined, the “HOLD” symbol stops flashing, and the weight is displayed on the screen for 2 minutes.

NOTE: If the weight cannot be automatically determined within 8 seconds, the scale will calculate an average weight over 4 seconds (during this time, the weight display (3) on the screen will flash).

- The measured weight remains stored in the scale, and the device will automatically switch off after 2 minutes. To turn it on again, press and hold the ON/OFF button for 3 seconds.

MANUAL HOLD (MECHANICAL HOLD):

If the baby is excessively restless or moving continuously, you can activate the manual hold function by pressing the HOLD button. After 5 seconds, the “HOLD” symbol will stop flashing, and the average weight will be displayed on the screen for 2 minutes.

DEACTIVATING THE AUTOMATIC HOLD FUNCTION (AUTO HOLD):

If you wish to use the scale without the automatic hold function, it can be deactivated.

- Once the display shows “0.000 kg,” you can deactivate the AUTO HOLD function by pressing and holding the HOLD button for 3 seconds.
- As soon as “0.000 kg” is displayed, the baby scale is ready for weighing.
- Carefully place the baby in the center of the weighing platform.
- To capture the baby’s weight, press the HOLD button. The recording of the weight is confirmed on the display when the “HOLD” symbol flashes for 3 seconds.
- The weight will remain displayed for 30 seconds, after which the scale will switch off automatically. Press the ON/OFF button to perform another measurement.
- To switch off the scale manually, press and hold the ON/OFF button for 3 seconds.
- To reactivate the automatic hold function, press and hold the HOLD button for 3 seconds.

NOTE: When the “HOLD + ON” symbol appears on the display, the automatic hold function is enabled.

7.3 TARE FUNCTION

The weight of any initial load (such as a weighing mat or cover) can be tared by pressing the TARE button. This ensures that during subsequent measurements, only the actual weight of the baby will be displayed.

Use this function when weighing the baby on a mat and while the baby is wearing a diaper and/or clothing.

- Switch on the scale. The scale is ready for weighing as soon as the display shows “0.000 kg.”
- Place the mat/diaper/clothing on the scale.
- Press the TARE button to activate the tare function. The weight reading will be reset to zero, and “0.000” will appear on the display. The scale is now ready for weighing.
- Carefully place the baby on the scale. The baby’s weight will be shown on the display.

If you placed a mat on the scale after switching it on and you wish to remove it before weighing, you can also use the tare function. When you press the TARE button, the “TARE” symbol appears on the display. To reset the scale to “0.000 kg,” remove the mat and press the TARE button again.

7.4 USING THE TAPE MEASURE

- Pull out the tape measure on the top of the baby scale if you want to measure the baby.
- In order to reel the tape measure back into the baby scale, press the round button next to the tape measure on the bottom of the baby scale.

8. CLEANING AND MAINTENANCE

- The accuracy of the measurements and service life of the device depend on its careful handling: The baby scale should be cleaned from time to time.
- This should be done using a damp cloth and, if necessary, a small amount of detergent.
- Do not submerge the baby scale in water. Do not rinse the baby scale under running water. Do not clean the baby scale in the dish washer.
- Never use harsh solvents or cleaning products!
- If you do not plan to use the baby scale for a long period of time, we recommend that you store it in the original packaging in a dry environment.
- During storage, do not place any objects on the baby scale.
- Store the baby scale away from children and pets.
- Remove the batteries from the baby scale.

9. TROUBLESHOOTING GUIDE

| DISPLAY | CAUSE | SOLUTION |
|-------------------------------|---|---|
| EEEE | The surface on which you have placed the scale is not stable and/or flat. | Place the scale on a flat stable surface without carpet. Retare the scale and repeat the measurement. |
| | The maximum scale weight of 20 kg has been exceeded. | Do not exceed the limit on the scale. |
| Incorrect weight is displayed | The surface on which you have placed the scale is not stable and/or flat. | Place the scale on a flat stable surface without carpet. Retare the scale and repeat the measurement. |
| | The baby kicks its legs and moves too much. | Try to calm the baby so you can get accurate measurements. |
| | Incorrect scale reset. Incorrect zero point of scale. | Turn the scale off, turn it back on and wait for "0.000 kg" to appear on the display. Repeat the measurement. |
| Lo | The batteries are empty. | Replace the batteries in the baby scale. |

WARNING!

If the device does not operate even after performing the above procedures, please contact the store where you purchased it or the nearest authorized distributor. Never turn on the device if any part of it is damaged or not functioning properly. The device should only be repaired by an authorized customer service representative or supplier.

By law, empty batteries must be disposed of at designated locations, recycling points or electronics retailers. The codes below are printed on batteries containing harmful substances:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury

For environmental reasons, do not dispose of the device in household waste at the end of its service life. Dispose of the device at a collection or recycling point.

This appliance complies with the following documents:
 Directive 1999/5/EC on radio equipment and telecommunications terminal equipment.
 Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment.

Manufactured for CANGAROO in PRC;

Manufacturer and Importer: Moni Trade LTD

Address: Bulgaria, Sofia, Trebich, Dolo str. 1 Tel.: +359 2 936 07 90 Website: www.cangaroo-bg.com

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